

# ROBOT MODELING AND CONTROL

ECE470S

## LAB2 : Forward and Inverse Kinematics for the KUKA Robotic Arm

### 1 Purpose

The purpose of this lab is to adapt the results of lab 1 for implementation in a real robot. The objective is to make the robot arm draw on paper.

### 2 Introduction

The KUKA robotic arm depicted in Figure 1 is an articulated manipulator with a three-degree-of-freedom wrist and a gripper. Highlighted in Figure 1 are the axes of rotation of joints 1 to 6. A pencil, not displayed

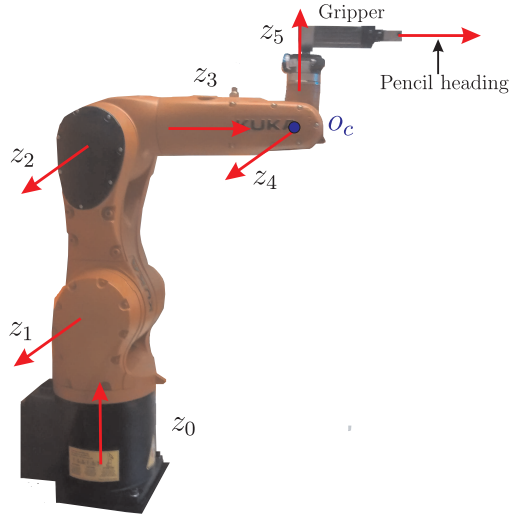


Figure 1: The KUKA robotic arm in the HOME configuration

in Figure 1, will be attached to the gripper and will point along the pencil heading direction indicated in the figure. The tip of the pencil is the end effector of the robot. The objective of this lab is to make the robot draw patterns on paper.

The kinematic parameters of the robot are detailed in Figure 2. The red elements in the figure are the joint axes. The values of various parameters are found in the table below.

$d_1$	400 mm
$a_1$	25 mm
$a_2$	315 mm
$a_3$	35 mm
$d_4$	365 mm
$a_6$	296.23 mm
$d_6$	161.44 mm

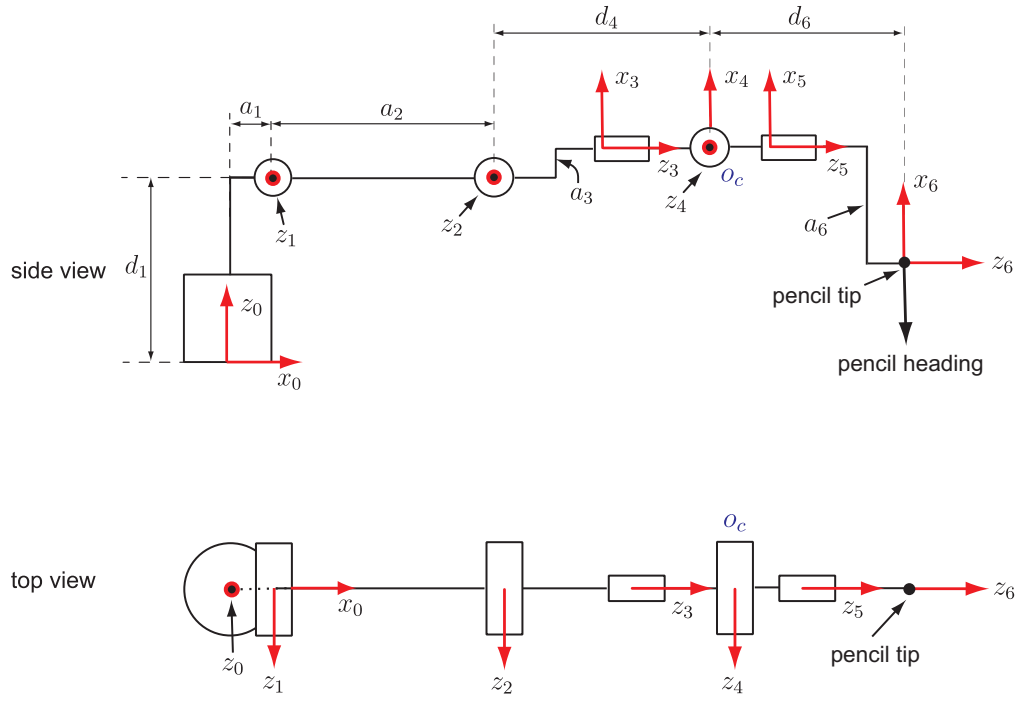


Figure 2: Schematic diagrams illustrating the KUKA kinematic parameters

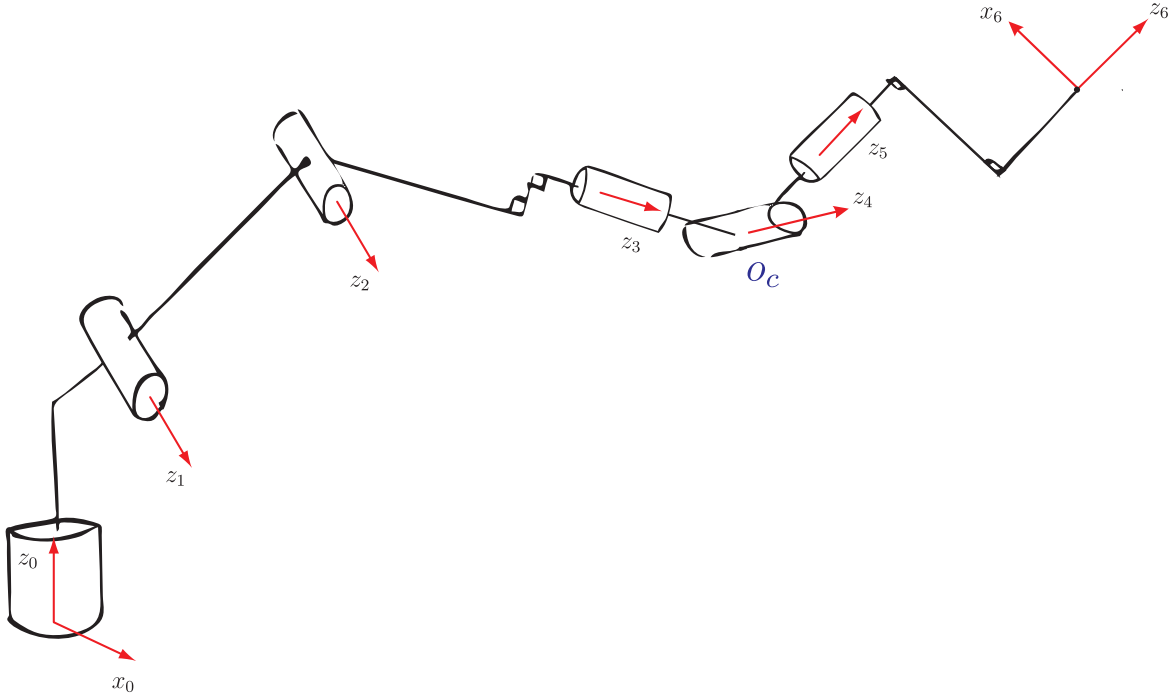


Figure 3: Schematic diagram of KUKA robot arm not at rest.

### 3 Preparation

Please submit a **complete** preparation at the beginning of the lab session. Prior to the lab, you should make sure that step 4 of the preparation returns consistent results. An incomplete or incorrect preparation will be penalized.

- Figure 3 depicts the robot in a generic configuration away from the rest position. Print out this page and draw on the figure coordinate frames  $o_1x_1y_1z_1$ ,  $o_2x_2y_2z_2$ ,  $o_3x_3y_3z_3$ ,  $o_4x_4y_4z_4$  and  $o_5x_5y_5z_5$  according to the DH convention, and write the DH table of the robot. In particular, choose the  $x_1$  axis parallel to the  $x_0$  axis,  $x_1 = x_0$  (this is done for convenience).

If your DH frame assignment is correct, you should find that when  $\theta_1 = \theta_2 = \theta_3 = 0$ , link 2 is parallel to the ground, and link 3 is vertical, pointing downward. A schematic diagram of the KUKA robot arm in the HOME configuration is illustrated in Figure 4 which you should verify corresponds to  $(\theta_1, \theta_2, \theta_3, \theta_4, \theta_5, \theta_6) = (0, \pi/2, 0, 0, \pi/2, 0)$ .

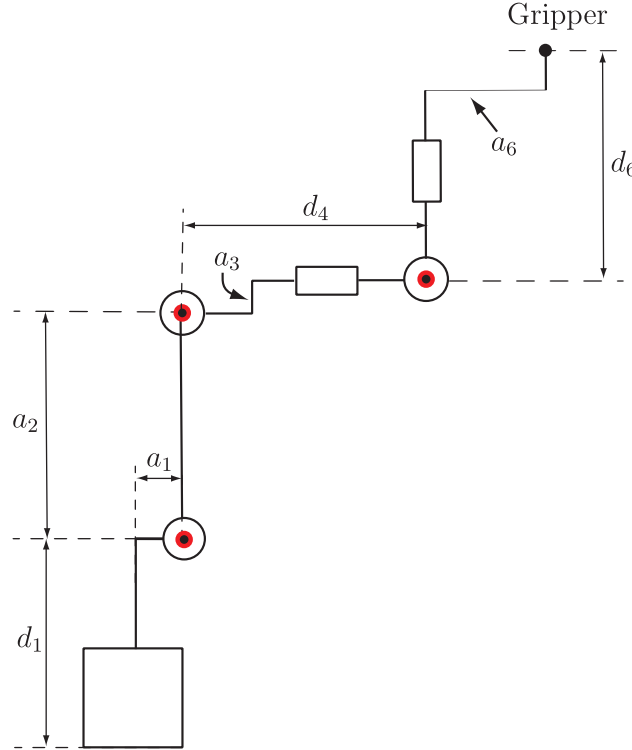


Figure 4: The KUKA robotic arm schematic diagram in the HOME configuration

- Derive the inverse kinematics of the robot. Specifically, given  $R_6^0(\theta_1, \dots, \theta_6) = R_d$  and  $o_6^0(\theta_1, \dots, \theta_6) = o_d^0$ , find  $(\theta_1, \theta_2, \theta_3, \theta_4, \theta_5, \theta_6)$ . You will find two solutions: elbow up and elbow down. Find the elbow up solution. Write your derivations neatly on paper. As in lab 1, one can solve the inverse kinematics problem by the technique of kinematic decoupling in which the problem is divided in two parts: inverse position and inverse orientation.

- The position of the wrist centre  $o_c$  is shown in Figure 2. First find  $(\theta_1, \theta_2, \theta_3)$  such that  $o_c^0(\theta_1, \theta_2, \theta_3) =$

$$\begin{bmatrix} x_c \\ y_c \\ z_c \end{bmatrix} = o_d^0 - R_d \begin{bmatrix} -a_6 \\ 0 \\ d_6 \end{bmatrix}.$$

- Then solve the equation

$$R_6^3(\theta_4, \theta_5, \theta_6) = (R_3^0)^\top R_d.$$

for  $(\theta_4, \theta_5, \theta_6)$ .

3. Modify your inverse kinematics function from lab 1 to incorporate the changes of this setup. Specifically, write Matlab functions `mykuka.m`, `forward_kuka.m`, and `inverse_kuka.m` as follows.

`myrobot = mykuka(DH)` defines the robot structure of the KUKA robot with the  $6 \times 4$  DH table you found earlier.

`H = forward_kuka(q,myrobot)` returns the homogeneous transformation matrix  $H$  of the end effector, where  $q$  is the  $6 \times 1$  vector of joint angles, and `myrobot` is the robot structure defined above.

`q = inverse_kuka(H,myrobot)` returns the  $6 \times 1$  vector of joint angles  $q = (\theta_1, \theta_2, \theta_3, \theta_4, \theta_5, \theta_6)$ , where  $H$  is the  $4 \times 4$  homogeneous transformation matrix  $H = \begin{bmatrix} R_d & o_d^0 \\ 0 & 1 \end{bmatrix}$ .

4. Test your software: you should get

```
>> kuka=mykuka(DH);  
>> forward_kuka([pi/5 pi/3 -pi/4 pi/4 pi/3 pi/4]',kuka)
```

ans =

0.1173	-0.3109	0.9432	368.9562
-0.8419	-0.5349	-0.0717	420.4832
0.5268	-0.7856	-0.3245	120.8570
0	0	0	1.0000

```
>> inverse_kuka(ans,kuka)
```

ans =

0.6283
1.0472
-0.7854
0.7854
1.0472
0.7854

## 4 Experiment

In this lab you will learn to interface the KUKA robot with Matlab, you will then calibrate the DH parameters by taking measurements, and finally you will use the Matlab functions you developed in your preparation to make the robot draw patterns on paper. We begin by recalling the Matlab commands that you can use to control the Kuka robot arm.

### 4.1 Safety procedures and manual operation of the robot

The KUKA robot you'll use in this course must be handled with great care, placing particular attention to issues of safety. Before you begin this lab, you need to have read Sections 1-4 of the document titled `Intro_KUKA_workstation`, and reviewed the `safety_checklist` document.

Confirm with your lab TA that you have understood the safety procedures and that you have read the document `Intro_KUKA_workstation`. Now perform the hands-on activity 1 in Section 4 of `Intro_KUKA_workstation`.

Before proceeding to the next section, show the TA what you have learned.

## 4.2 Review: Commanding the KUKA robot arm through Matlab

Below are the Matlab commands you will use to control the Kuka robot arm.

- `startConnection` establishes connection between Matlab and KUKA.
- `stopConnection` terminates connection between Matlab and KUKA.
- `getAngles()` returns the vector of current joint angles  $q = (\theta_1, \theta_2, \theta_3, \theta_4, \theta_5, \theta_6)$ .
- `stop()` stops KUKA motion.
- `moveAxis(axis,vel)` moves a single KUKA axis. `axis` takes a value between 1 and 6 that corresponds to the joint angle to be commanded. `vel` is a signed value that determines the commanded angular speed of the joint. In this lab `vel` should be set to no greater than 0.01.
- `setAngles(q,vel)` sets KUKA arm angles to those defined in `q` to within a small tolerance. `vel` corresponds to the speed of motion. For this lab, set `vel` to 0.04. To cancel the command during execution press `ctrl c` and run `stop()` in the Matlab command window.
- `setHome(vel)` sets KUKA arm to the HOME configuration in Figure 1. `vel` corresponds to the speed of motion. For this lab, set `vel` to 0.04. To cancel the command during execution press `ctrl c` and run `stop()` in the Matlab command window.
- `setGripper(state)` sets the gripper state. `setGripper(0)` closes the gripper; `setgripper(1)` opens the gripper.

**Warning:** Never use the `clear all` command in Matlab. Instead use `clearvars -except udpObj` to avoid errors from occurring.

The steps to **connect KUKA to the external PC**:

- Run `startConnection` in Matlab.
- On the KUKA SmartPad run `RSI.Ethernet.src` until the line `RSI.MOVECORR()`.
- To run the program, hold down half-way one of the enable buttons on the back of the KUKA SmartPAD. The enable buttons are labelled 3 and 5 in Figure 5. If you press too hard, an emergency stop will be activated. While holding down the enable button, press and hold the run button labelled 10 in Figure 6. Note that when running the line `RSI.ON` a warning will appear saying ‘Caution - sensor correction is activated’. Simply confirm this warning to continue.

The steps to **command the KUKA arm from Matlab** are as follows:

- Run the desired command in the Matlab command window. The robot will not yet move.
- Hold down half-way one of the enable buttons on the back of the KUKA SmartPAD.
- While holding down the enable button, press and hold the run button. Now the robot will carry out the desired command as long as both the enable and run buttons are being held down.
- To make the robot stop moving, simply release the enable button. To cancel the current command, press `ctrl c` and run `stop()`.
- If the robot ever hits an object while moving (such as the ground), immediately release the enable button and the robot will stop. Then press `ctrl c` in the Matlab control window to cancel the command and run `stop()`. Immediately call a TA to resolve the collision.

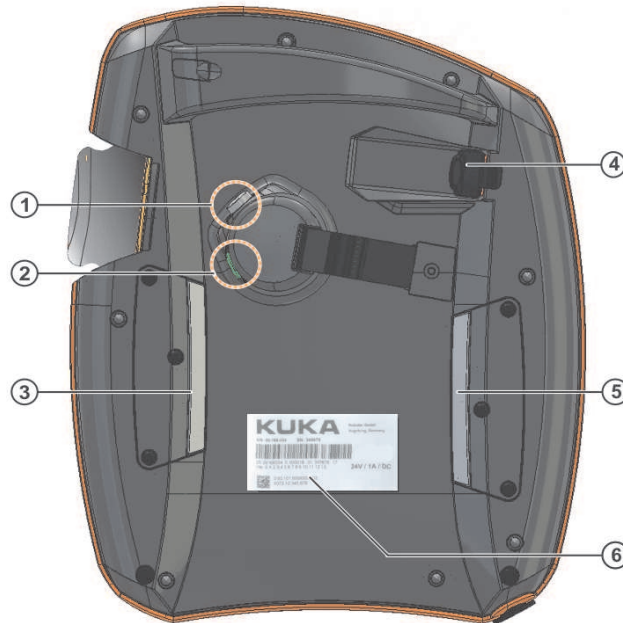


Figure 5: 3 and 5 are Enable buttons on the KUKA SmartPAD.

- If the robot does not move when commanded, it means an error has occurred. Call a TA to resolve the issue.
- For safety, students must never open the safety gate enclosing the robot unless instructed to do so. The enable and start button must never be pressed while the safety gate is open.

**Warning:** Make sure to **FIRST** connect KUKA to Matlab using `startConnection.m`. **THEN** run the RSI from the Kuka SmartPad

The steps to **disconnect KUKA from the external PC**:

- When finished running the lab, return to the program `RSI.Ethernet.src` on the SmartPAD, click on line 32 and then click **Block selection** to set the program cursor to the command `ret=RSI_OFF()`. Then continue running the program to the end. The program status indicator will turn black to indicate the program is complete.
- Run `stopConnection` in Matlab.

The steps to **deal with an RSI connection problem**:

- Type `CTRL+C` within the Matlab environment.
- There are two options available on the top portion of the SmartPAD screen: **cancel** or **reset**. Use **reset**, not cancel, to reinitialize the RSI program from the Kuka SmartPAD.
- After resetting the RSI program, push the play button multiple times until the program reaches the **RSI\_MOVECORR()** line.
- In Matlab, run `getAngles()` to verify that the RSI connection has been correctly re-established.

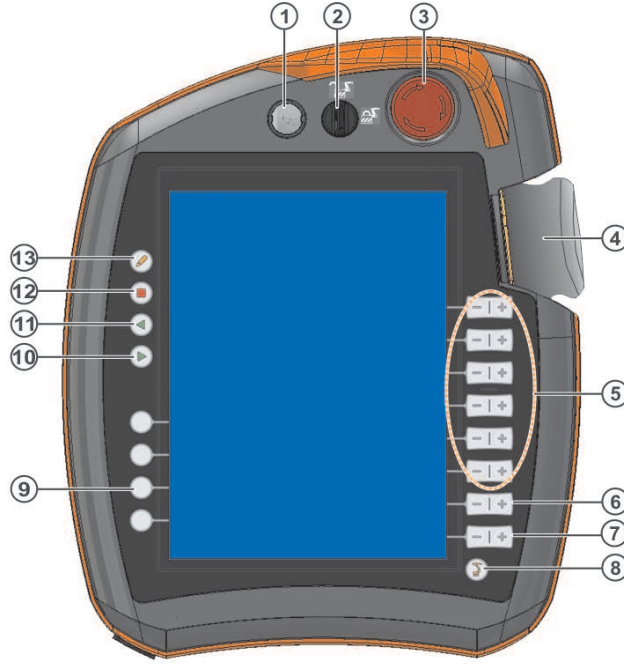


Figure 6: 10 is the Run button on the KUKA SmartPAD.

### 4.3 Calibration of DH parameters

The DH parameters provided to you in the preparation were approximate. You need to calibrate them in order to improve the accuracy of the forward kinematics function. The main source of the inaccuracy are the parameters  $a_6$ , and  $d_6$ . We will tune these using a simple optimization approach. These are the steps we will perform to calibrate the robot:

- **Data collection.** We will collect data pairs  $(Q_i, X_i)$  for three sample points, where  $Q_i$  is the joint vector of the robot, and  $X_i$  is the corresponding position of the pencil tip in frame 0.
- **Parametric robot structure.** We will parametrize the robot structure definition in `mykuka.m` by two parameters,  $\delta_1$  and  $\delta_2$ , which will perturb the DH parameters  $a_6$  and  $d_6$ . The parameters  $\delta_1$  and  $\delta_2$  are placed in a  $2 \times 1$  vector `delta`.
- **Cost function.** We will define a cost function `deltajoint.m` (provided to you) that, given the vector `delta` of DH parameter perturbations, uses forward kinematics to predict coordinates  $\hat{X}_i$  of the end effector corresponding to the joint vectors  $Q_i$  measured in step 1. The cost measures the discrepancy between the predictions  $\hat{X}_i$  and the actual  $X_i$  measured in step 1.
- **Cost minimization.** Using the Matlab command `fminunc.m`, we will find the vector `delta` minimizing the above cost function. This vector constitutes the perturbations to the parameters  $a_6$  and  $d_6$  that minimize the discrepancy between the measured  $X_i$  and the  $\hat{X}_i$  predicted using forward kinematics.
- **Calibrated robot structure.** Using `delta` found in the cost minimization step, we will define a new, calibrated robot structure.

Now the detailed steps you need to perform.

1. Place the robot's end effector (pencil) at three different sample points on the floor using the jog keys of the KUKA SmartPad. *Mark one of the three sample points on the grid paper.* We will need this marked point later.

Make sure the points are not too close to each other. For each sample point, perform the following measurements:

- Use the command `getAngles()` to read the joint angles of the robot at the sample point, and store the returned value in a vector  $Q_i$ , where  $i = 1, 2, 3$  is the index of the sample point.
- On the SmartPAD, read the  $(x, y, z)$  coordinates of the tool tip in frame 0, and save them in a vector  $X_i$ ,  $i = 1, 2, 3$ . Make sure to choose “D\_pen” as your reference frame on the KUKA Smart-Pad. See the appendix for instructions on how to do this. Make sure to note these coordinates in the same units as those used in your DH table (millimetres).  
You have thus obtained these measurements:  $Q_i := [\theta_1^i, \theta_2^i, \theta_3^i, \theta_4^i, \theta_5^i, \theta_6^i]^\top$  and  $X_i := [x_i, y_i, z_i]^\top$ , for  $i = 1, 2, 3$ .
- Using the vectors just collected, update variables `X1`, `X2`, `X3`, `Q1`, `Q2`, `Q3` in the `deltajoint.m` Matlab file (lines 6-11).
- Update variables `X1`, `X2`, `X3` in the `FrameTransformation.m` Matlab file (lines 4-6). You will use this function in the next two sections.

2. Duplicate your `mykuka.m` function and save it as `mykuka_search.m`. This new function should work like this:

```
myrobot = mykuka_search(delta)
```

where  $\delta$  is a  $2 \times 1$  vector containing the perturbations of the parameters  $a_6$ , and  $d_6$ . In other words, the DH table in `mykuka_search.m` should contain parameters  $a_6 + \delta(1)$ , and  $d_6 + \delta(2)$ .

3. Open and inspect the provided Matlab function `deltajoint.m`:

```
X_error = deltajoint(delta)
```

This function uses `mykuka_search.m` and `forward_kuka.m` to compute the sum of the errors between the measured vectors  $X_i$  and their estimates computed using forward kinematics. Specifically,

```
H1=forward_kuka(Q1,kuka);
H2=forward_kuka(Q2,kuka);
H3=forward_kuka(Q3,kuka);
```

```
X_error=norm(H1(1:3,4)-X1)+norm(H2(1:3,4)-X2)+norm(H3(1:3,4)-X3);
```

You don't have to make any change to this function, but make sure you understand its operation.

4. Use the Matlab command `delta = fminunc(@deltajoint,[0 0])` to find the optimal vector `delta` which minimizes the total joint variable error computed through the function `deltajoint.m`.
5. Redefine the robot structure using the updated DH parameters. You can do this by issuing the Matlab command `myrobot = mykuka_search(delta)`, where `delta` is the vector of parameters you obtained at the previous step.
6. Test your calibration. Command the robot to the HOME configuration. Solve an inverse kinematics problem `q = inverse_kuka(H,myrobot)`, where the rotation matrix part of `H` is

$$R_6^0 = \begin{bmatrix} 0 & 0 & 1 \\ 0 & -1 & 0 \\ 1 & 0 & 0 \end{bmatrix}$$



and the translation part is  $o_6^9 = \mathbf{X1}$ , where  $\mathbf{X1}$  is the sample point from step 1 that you have marked on the grid paper. Then using `setAngles(q,0.04)`, verify that the effector goes to the test point on the grid paper, and that the pencil is pointing vertically downwards.

Note down the calibrated DH parameters  $a_6$  and  $d_6$  that you found. You will reuse them in Lab 4.

**Show these results to your lab TA.**

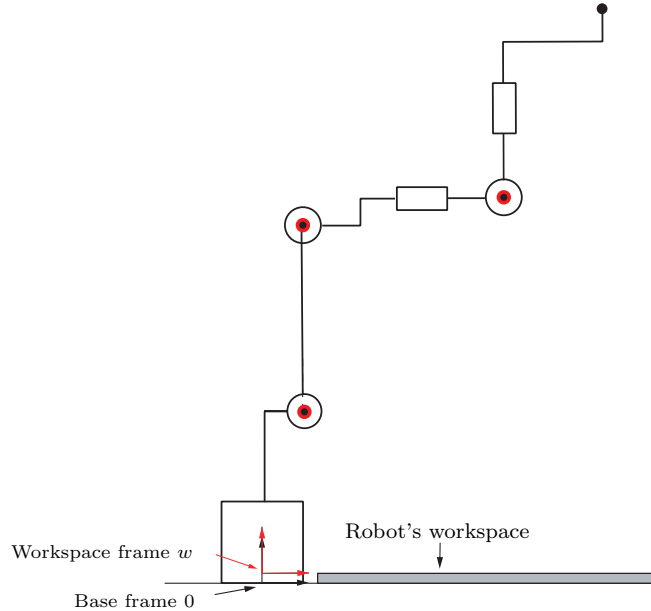


Figure 7: Base and Workspace frames of the Kuka robot.

#### 4.4 Workspace frame versus Base frame

The base frame of the Kuka robot, displayed in Figure 7, is assigned by the manufacturer. The workspace of the robot setup at the University of Toronto is the surface of a wood platform. This platform has a certain thickness and is not perfectly horizontal. For this reason, we make use of a *Workspace frame* which we label frame  $w$ , displayed in Figure 7, whose  $z$  axis is close to being vertical, but not perfectly so, and whose origin is higher than the origin of the Base frame to reflect the thickness of the Workspace platform.

The Matlab function `FrameTransformation.m` provided to you utilizes the  $(Q_i, X_i)$  pairs you found in the previous section to determine the homogeneous transformation matrix  $H_w^0$  that converts points from the coordinates of frame  $w$  to those of frame 0.

In Section 4.5 you will draw patterns on a sheet of paper taped on the robot's Workspace. The specifications that will be given to you will be expressed in the coordinates of frame  $w$ , and you will use `FrameTransformation.m` to convert them to frame 0. To get ready for drawing patterns, we now practice the conversion in question. To this end, we want to make the robot's pencil reach the point  $p^w = [600 \ 100 \ 10]^T$  (units are in millimetres). Before going forward, make sure to update variables  $\mathbf{X1}$ ,  $\mathbf{X2}$ ,  $\mathbf{X3}$  in the `FrameTransformation.m` Matlab file (lines 4-6), as detailed in Step 1 of Section 4.3.

1. Define the target point in Workspace coordinates: `p_workspace = [600; 100; 10]`.
2. Convert the point in base frame coordinates: `p_baseframe = FrameTransformation(p_workspace)`.
3. Choose a desired orientation of the end effector with respect to the base frame (whereby the pencil points vertically downward): `R = [0 0 1; 0 -1 0; 1 0 0]`.

4. Define the desired homogeneous transformation matrix of the end effector with respect to frame 0:  $H = [R \ p\_baseframe; \ zeros(1,3) \ 1]$ .
5. Use inverse kinematic to find the desired joint variables:  $q = \text{inverse\_kuka}(H, \text{myrobot})$ .
6. Command the robot to the target position:  $\text{setangles}(q, 0.04)$ .

## 4.5 Draw Patterns

Now that you have calibrated the robot, you are ready to draw patterns. The idea is to generate desired end effector trajectories, use the `inverse_kuka.m` function to translate them into desired joint trajectories, and then commanding these joint trajectories to the robot. Throughout this section, we will want to keep the pen vertical, pointing downward. This translates into requiring that  $R_6^0$  be given as

$$R_6^0 = \begin{bmatrix} 0 & 0 & 1 \\ 0 & -1 & 0 \\ 1 & 0 & 0 \end{bmatrix}. \quad (1)$$

In the development that follows, you will make  $o_6^0$  follow paths corresponding to a number of patterns. Meanwhile,  $R_6^0$  will be the constant constant matrix in (1).

1. Write a Matlab function `mysegment.m` performing these tasks. Generate a  $3 \times 100$  matrix `X_workspace` containing 100 end effector positions in the coordinates of frame  $w$  describing a straight line segment on the table parallel to the  $y^0$  axis. Specifically, for  $i = 1, \dots, 100$ , `X_workspace(:, i)` is a  $3 \times 1$  vector of the form  $[\bar{x}_i \ \bar{y}_i \ \bar{z}_i]^T$ , where  $\bar{x}_i = 620$  mm,  $\bar{y}_i$  ranges from  $-100$  mm to  $100$  mm and  $\bar{z}_i = -1$  mm (this is to guarantee that the pencil makes contact with the paper).

Generate a  $3 \times 100$  matrix `X_baseframe` whose columns are the columns of `X_workspace` converted to the coordinates of frame 0 using the function `FrameTransformation.m`, as described in Section 4.4.

Next, for each column `X_baseframe(:, i)` of `X_baseframe`, form a matrix  $H$  using  $R_6^0$  as in (1), and  $o_6^0$  equal to `X_baseframe(:, i)`. Compute the joint angles using `inverse_kuka.m` and command them to the robot using `setangles` in a loop.

Note: you may need to slightly adjust the  $\bar{z}_i$  coordinates in the matrix `X_workspace` to guarantee that the pencil pushes down hard enough on the surface to create a visible line.

**Show your results to the lab TA. Take a photo of the resulting pattern drawn by the robot.**

2. Write a Matlab function `mycircle.m` analogous to `mysegment.m`, but now such that the end effector draws a circle of radius 50 mm centred at the point  $[620 \ 0 \ -1]^T$  in the workspace frame.

**Show your results to the lab TA. Take a photo of the resulting pattern drawn by the robot).**

3. Invent a creative pattern, write a Matlab function for it, and make the robot draw it. Take a photo of the result, and shoot a movie of the robot in action. For instance, try executing the pattern `jug.xlsx` posted on Blackboard scaled by ten and centred about the origin of the frame marked on the grid paper. Begin with the commands:

```
data=xlsread('jug.xlsx');
xdata=550 + 10*data(:,1);
ydata=10*data(:,2);
zdata=-ones(length(data),1);
```

Can you make something better?


## 5 Submission

Each student will submit on Quercus a preparation report *in advance of the lab*. No later than one week after your lab session, submit on Quercus a zipped folder containing:


1. Your Matlab functions executing the steps in Section 4. The main file should be called **Lab2.m**. The lab code should be thoroughly commented.
2. Photos of the three patterns drawn by the robot (segment, circle, and your pattern).
3. One movie of the robot drawing the creative pattern you came up with.

### Appendix: How to read the coordinates of the pencil tip on the SmartPAD


1. Open the main menu ((A) in the figure below) on the SmartPAD, and select Display > Actual position ((B) in the figure below).

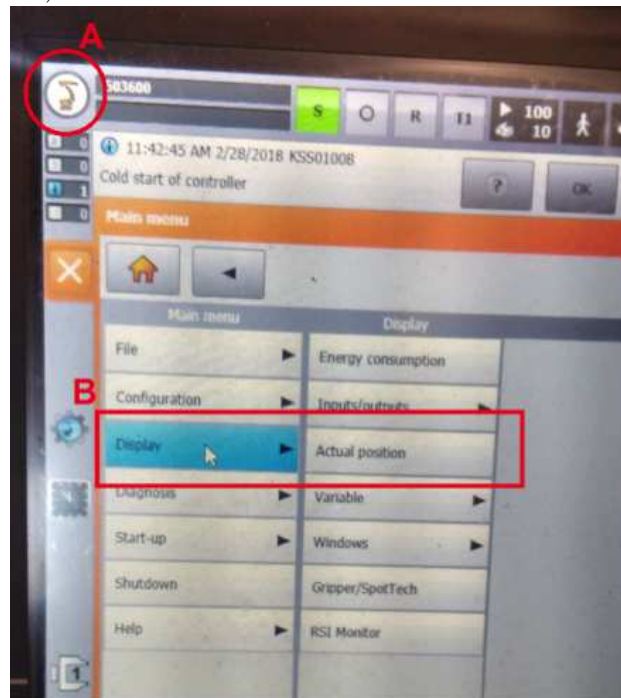

2. You should see a display showing the coordinates ( (C) in the figure below). The “Cartesian Position” option should be selected (this can be toggled with “Axis Specific”, see (D)). Verify that the tool selected is **pencil** and that the frame is **\$NULLFRAME** (see (E)). Currently, on all machines this is set to tool 11 and base 0. *If the incorrect tool is selected, the coordinates will not be accurate.*

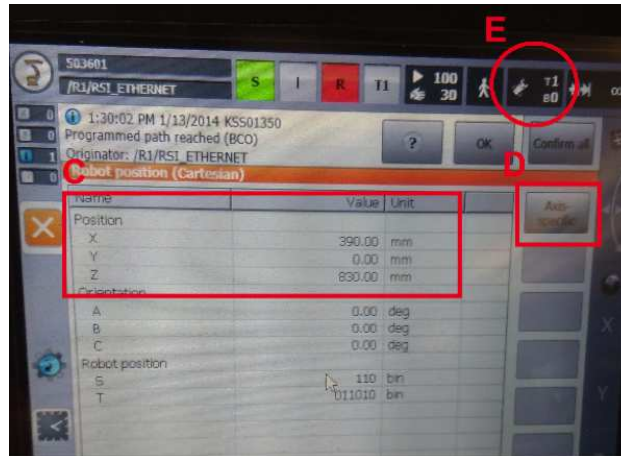

- The tool/base should be as shown in the image below. Note that the “Flange” option should be selected.

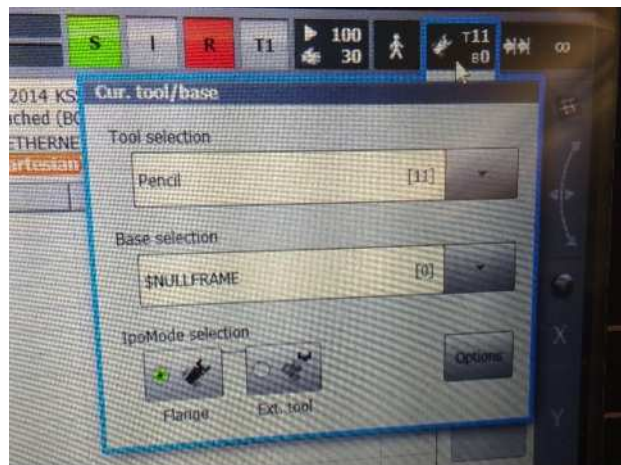

- If `RSI_ETHERNET.src` is cancelled/rerun, the tool always defaults to tool [1] . After running `RSI_ETHERNET.src`, always change the tool back to the pencil tool.